SciCan Data *Logger*

- **Operator's Manual**
- **Manuel de l'opérateur**
- **Benutzerhandbuch**
- **Manuale per l'operatore**
- **Manual del operador**

SciCan Data Logger

Installing the SciCan Data Logger on your Statim	8
Data Logger Troubleshooting for Statim	12
Installing the Data Logger on your Hydrim	13
Data Logger Troubleshooting for Hydrim	14
Specifications	15

Statim and Hydrim are registered trademarks of SciCan Ltd. STAT-DRI is a trademark of SciCan Ltd.
All other trademarks referred to in this manual are the property of their respective owners.

For all service and repair inquiries:

In Canada	1-800-870-7777
United States:	1-800-572-1211
Europe:	49-0821-278-930
International:	(416) 446-4500
Email:	techservice.ca@scican.com

Manufactured by: SciCan Ltd. 1440 Don Mills Road, Toronto ON M3B 3P9 CANADA Phone: (416) 445-1600 Fax: (416) 445-2727 Toll free: 1-800-667-7733 	SciCan Inc. 701 Technology Drive Canonsburg, PA 15317 USA Phone: +1 724 820 1600 Toll free: +1 724 820 1479 Fax: 1-888-632-6031
--	---

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SciCan Ltd.

1440 Don Mills Road,
Toronto ON M3B 3P9
CANADA

Teléfono: (416) 445-1600

Fax: (416) 445-2727

Llamada gratis: 1-800-667-7733

SciCan Inc.

701 Technology Drive
Canonsburg, PA 15317
EE.UU.

Teléfono: +1 724 820 1600

Llamada gratis: +1 724 820 1479

Fax: 1-888-632-6031

Representante de la UE:

SciCan GmbH

Wangener Strasse 78

88299 Leutkirch

ALEMANIA

Teléfono: +49 (0)7561 98343 - 0

Fax: +49 (0)7561 98343 - 699

SciCan Medtech AG

Alpenstrasse 16, 6300 Zug Suiza

Teléfono: (41-41) 727-70-27

Fax: (41-41) 727-70-29


INSTALLING THE SCICAN DATA LOGGER ON YOUR STATIM

SciCan's Data Logger can record and store cycle information onto a mass storage device (MSD) such as a USB Flash Drive or SD memory card.

These instructions are for units with a 9-pin communication port (RS232) located on the back of the unit.

Follow the steps below before connecting the Data Logger.

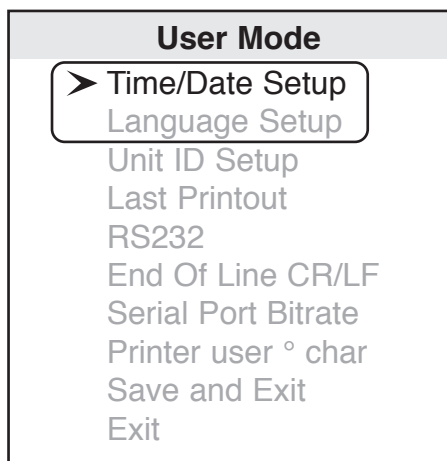
User Setup Menu

To install the Data Logger, access the user setup menu by powering up the unit while pressing the **STOP** button. Follow the instructions below to complete installation steps.

Step 1

Selecting the USB Flash/MSD Option

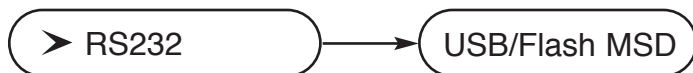1. Power up unit while pressing the **STOP** button to access the User Setup menu.
2. Using the UNWRAPPED and WRAPPED buttons, scroll to RS232 and select it by pressing the RUBBER AND PLASTIC button.
3. From the RS232 menu, use the UNWRAPPED and WRAPPED buttons to scroll down to the USB/Flash MSD option and press the RUBBER AND PLASTIC button to select and return to the User Setup menu.

Setting the Serial Port Bitrte to 9600

➤ Serial Port Bitrte

9600

1. From the User Setup menu, use the UNWRAPPED and WRAPPED buttons to scroll to the Serial Port Bitrte menu and select it using the RUBBER AND PLASTIC button.
2. From the Serial Port Bitrte menu, use the UNWRAPPED and WRAPPED buttons to scroll down to 9600 and press the RUBBER AND PLASTIC button to select and return to the User Setup menu.

Setting the Printer user ° Character; (example: 134°C)

➤ Printer user ° char

32 [0x20]

1. From the User Setup menu, use the UNWRAPPED and WRAPPED buttons to scroll to the Printer user ° char menu and select it using the RUBBER AND PLASTIC button.
2. From the Printer user ° Char menu, using the UNWRAPPED button to increase the displayed value by one, and the WRAPPED button to increase that value by ten, enter the value 32 [0x20] and press the RUBBER AND PLASTIC button to accept and return to the User Setup menu.

Save and Exit

➤ Save and Exit

A Save and Exit must be performed after the above settings are completed. If this is not done, the information will revert to its prior settings.

1. From the User Setup menu, use the UNWRAPPED and WRAPPED buttons to scroll to Save and Exit.
2. Select it by pressing the RUBBER AND PLASTIC button.


The LCD will display time and date information and a sequence of messages:

HH:MM DD/MM/YYYY

“MSD NOT DETECTED”/“INSERT
MSD/FLASH”/“SELECT A CYCLE”

Step 2

Selecting the Time and Date

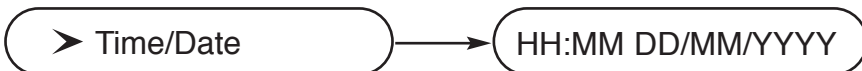1. Power up unit while pressing the **STOP** button to access the User Setup menu.
2. Using UNWRAPPED and WRAPPED buttons, scroll to Time/Date Setup and press the RUBBER AND PLASTIC button to select it.
3. From the Time/Date Setup menu, set the time and date using the UNWRAPPED and WRAPPED buttons to change the displayed values and the RUBBER AND PLASTIC button to select. Press the **STOP** button when completed.

Step 3

Setting Unit ID

NOTE: If the unit ID was set according to the Statim Operator's Manual, this procedure need not be repeated.

1. Power up unit while pressing the **STOP** button to access the User Setup menu.
2. Using the UNWRAPPED and WRAPPED buttons, scroll to the Unit ID Setup and press the RUBBER AND PLASTIC button to select it.
3. From the Setup Unit ID menu, use the UNWRAPPED and WRAPPED buttons to change the displayed values and the RUBBER AND PLASTIC button to select and move to the next digit. Press the **STOP** button when completed.

Step 4

Connecting the SciCan Data Logger

1. Ensure that both the Statim unit and the SciCan Data Logger are off.
2. Connect the SciCan Data Logger to the Statim unit by using the serial cable.

3. Power up the SciCan Data Logger
4. Power up the Statim unit.
5. The LCD will display the following sequence of messages:

HH:MM DD/MM/YYYY
“MSD NOT DETECTED”/”INSERT
MSD/FLASH”/”SELECT A CYCLE”

6. Insert USB flash drive or SD memory card.
7. After a few seconds the LCD will display the following sequence of messages:

HH:MM DD/MM/YYYY
USB/FLASH DETECTED/SAFELY
REMOVE MSD/SELECT A CYCLE

DATA LOGGER TROUBLESHOOTING FOR STATIM

If display reads: MSD NOT CONNECTED <> INSERT MSD/FLASH

Follow the steps below to ensure everything is setup properly.

1. Check the serial cable connection.
2. Check the power connection.
3. Ensure the lower red LED is lit. It is located between the Data Logger's serial port and power input.
4. Check that the Mass Storage Device (MSD) is properly inserted in the appropriate slot of the Data Logger.
5. Repeat the instructions for **Installing the Data Logger on your Statim.**

If display reads: SAFELY REMOVE MSD <> MSD/FLASH DETECTED

The MSD/Flash can be safely removed without affecting the data.

If display reads: MSD/FLASH FULL <> REPLACE MSD

The MSD is full. Export the data.

If MSD has: Missing lines of Data

Refer to the STATIM screen to confirm successful sterilization.

Reset the Data Logger by unplugging its power, disconnecting the MSD and waiting 10 seconds. Then re-connect the power adapter and insert the MSD into the Data Logger.

If the problem persists, contact the SciCan Service Center.

If MSD has: Corrupt or unreadable files and/or directories

Refer to the STATIM screen to confirm successful sterilization.

The MSD may have been unplugged while data was being written to it. The MSD should not be unplugged until after "SAFELY REMOVE MSD <> MSD/FLASH DETECTED" is displayed.

The corrupted files or directories may be lost. Reformat the MSD on your computer.

INSTALLING THE DATA LOGGER ON YOUR HYDRIM

The Hydrim C51wd has a 9-pin communications port (RS232) the back, and can be used with an external printer or the Data Logger. To use the Data Logger, follow these steps:

Turn off the Hydrim and Data Logger before connecting the Data Logger.

1. With the Data Logger connected and powered on, turn on the Hydrim and press the *i* on the bottom of the touchscreen to move to the Menu screen.
2. In the Setup menu, select Printer Selection.
Optional step. Setting Date and Time if not previously set. In the Setup menu select Date/Time and then Date Setup. Enter current date. When finished, press back arrow to return to Date/Time menu Select Time Setup. Enter current time. When finished, press back arrow to return to Setup menu.
3. Select USB Flash/MSD. Press the back arrow to return to the Setup Menu.
4. In the Setup Menu, select Baud rate.
5. Select 9600 as the Baud rate. Use the back arrow to return to the Start screen.
6. Now the Hydrim will write its cycle information to the memory stick inserted in the Data Logger.

DATA LOGGER TROUBLESHOOTING FOR HYDRIM

If display reads: “PRINTER FAULT”

Follow the steps below to ensure everything is setup properly.

1. Check the serial cable connection.
2. Check the power connection.
3. Ensure the lower red LED is lit. It is located between the Data Logger’s serial port and power input.
4. Check that the Mass Storage Device (MSD) is properly inserted in the appropriate slot of the Data Logger.
5. Repeat the instructions for **Installing the Data Logger on your Hydrim.**
6. Check to MSD to ensure it is not full. If MSD is full, export the data.

If display reads: Missing lines of Data on MSD/FLASH

Refer to the Hydrim screen to confirm successful sterilization.

Reset the Data Logger by unplugging its power, disconnecting the MSD and waiting 10 seconds. Then re-connect the power adapter and insert the MSD into the Data Logger.

If the problem persists, contact the SciCan Service Center.

If display reads: The file or directory in the MSD is corrupted or unreadable.

Refer to the Hydrim screen to confirm successful sterilization.

The MSD may have been unplugged while data was being written to it. The MSD should not be unplugged until after “SAFELY REMOVE MSD <> MSD/FLASH DETECTED” is displayed.

The corrupted files or directories may be lost. Reformat the MSD on your computer.

SPECIFICATIONS

Connection:	Serial Cable DB9F-DB25M
Baud Rate:	9600
Electrical Rating:	100V-240V, 0.6A, 50/60 Hz, multiplug
Env. Operating Temperature:	0-40°C
Env. Operating Relative Humidity:	20-85% (non-condensing)

ACKNOWLEDGED MSD

USB Flash Drives:

- Cruzer mini 1 GB
- Memorex Travel drive 1 GB
- Kingston Travel drive 512 MB
- PNY Attache 1 GB
- Lexar Jump Drive 1 GB
- SanDisk 2 GB

SD Memory Cards:

- Kingston 512 MB
- Lexar 512 MB
- PNY 512 MB
- Ultra 512 MB
- Kingston Ultimate 512 MB

Other USB Flash Drives and SD Cards may also be compatible with this Data Logger.

Recommendations:

SciCan has testing and recommends using the following printers as an external printer option, that can be purchased from your local computer store:

Suggested External Printers by SciCan	End Of Line CR/LF	Serial Port Bitrate	Printer user ° char
Epson TM-U220D (C31C515603)	CR/LF	9600	248 [0xF8]
Citizen IDP-3110-40 RF 120B	CR	9600	N/A
Star Micro SP212FD42-120	CR	9600	210 [0xd2]
Star Micro SP216FD41-120	CR/LF	9600	210 [0xd2]
Star Micro SP512MD42-R	CR/LF	9600	210 [0xd2]

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